

FINAL REPORT
(DURIP 00) Scanning Electron Microscope (SEM)

AFOSR Award# F49620-00-1-0261
UD Account# 3-3-21-2570-03

April 1, 2000 – March 31, 2001

George C. Hadjipanayis
Principal Investigator


Submitted to:
Dr. Harold Weinstock
AFOSR/NE
801 N. Randolph Street
Room 732
Arlington, VA 22203-1977
(703) 696-8572

20011003 114


The purchase of a high power field emission Scanning Electron Microscope under DURIP FY2000 grant F49620-00-1-0261 has enhanced the Magnetism lab microanalysis capabilities. The Department of Physics at the University of Delaware purchased a state-of-the-art JEOL JSM 6335F Scanning Electron Microscope in the fall of 2000. The purchase price of \$235,000.00 was negotiated to include a high resolution Backscatter detector, motorized X and Y stage control, and infrared specimen chamber scope. The installation was completed in a timely manner and the instrument was turned over to the Magnetism Lab in November 2000.

The increased analytical capabilities of the Magnetism Lab, since the purchase of the field emission SEM, has made the observation of nanostructures in hard and soft magnetic material routine. Having a high power field emission SEM means the laborious sample preparation and extensive time required for TEM imaging is no longer required. Typically, the resolution of the SEM (1.5nm at 15KV) and its depth of field and microanalysis capabilities provide detailed information of samples that previously required sending them to various other facilities. One of the unique features of the field emission SEM is the ability to view elemental maps using EDS techniques. The relatively small beam size provides the resolution to map microstructures and to identify inhomogeneities in samples in the range of a few nanometers. Since the installation of the JSM 6335F, our research activities have increased. Quick screening of samples and their associated microstructure can now be performed in a timely manner.

Actual JSM6335F images

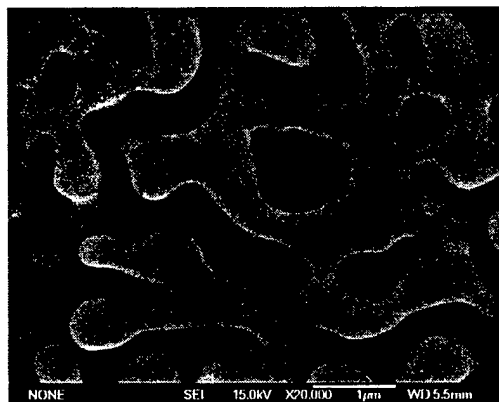

Porous Gold Microstructure

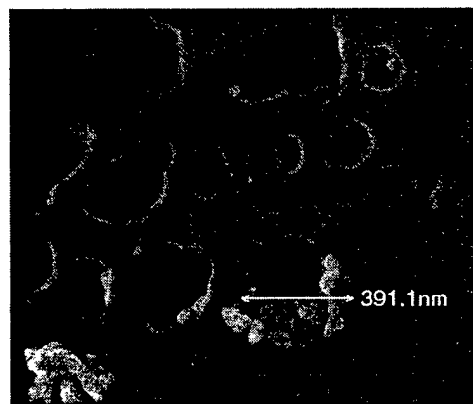

Silver Nanoparticles in
Polystyrene

One of the more recent undertakings of the UD magnetism lab has been the production of nanoparticles in a matrix using a novel sputtering based deposition technique. Shown below are SEM images of Fe nanoparticles in a BN matrix. Images at 30k magnification give particle sizes in the range of 20~100nm. At higher magnification (330k), it was possible to investigate clustering of nanoparticles. The image clearly shows the separation of the nanoparticles of a cluster enveloped by a possible passivation layer.

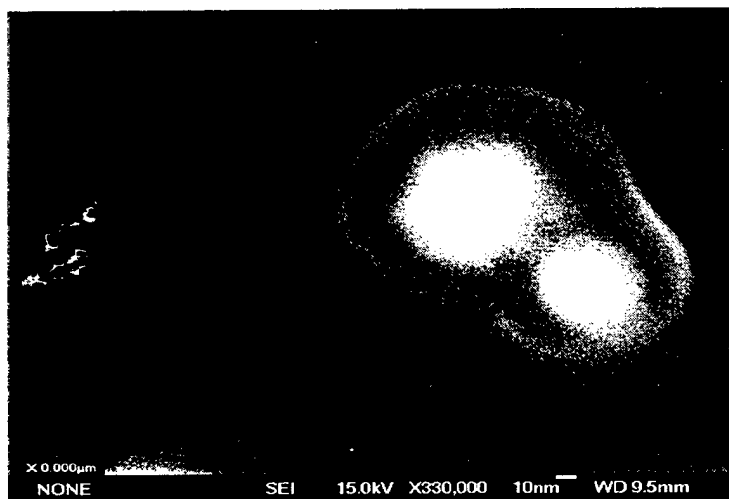
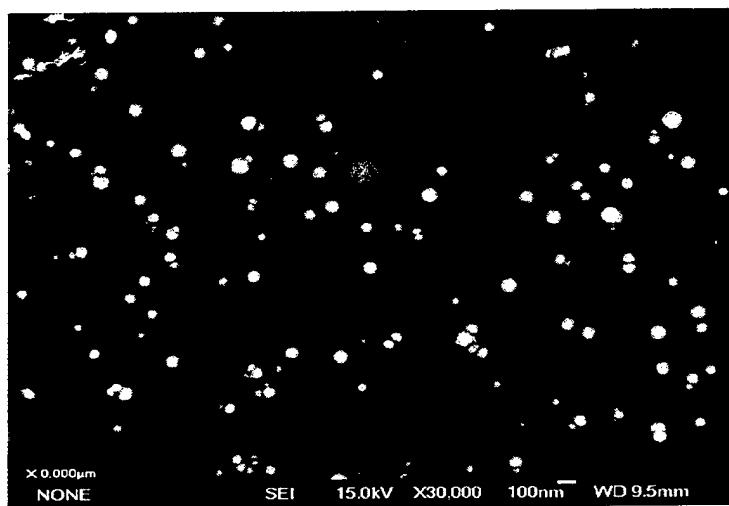


Figure 1 shows a typical SEM microstructure of cast SmCo_2Cu_3 magnets with EDS chemical composition analysis. There are many black dots in the right grain in figure 1a. EDS results found these dots had the same composition as matrix-1:5 phase, which means that these dots are originally from matrix phase. In higher magnification SEM images, some irregular particles were observed near the grain boundary (figure 1b). From EDS chemical composition analysis (table1), particle 1, particle 5, and particle 6 are the pieces of matrix phase, whereas particle 2, particle 3, and particle 4 show very high Sm content and considerable O content, which indicates that these particles are Sm oxides.

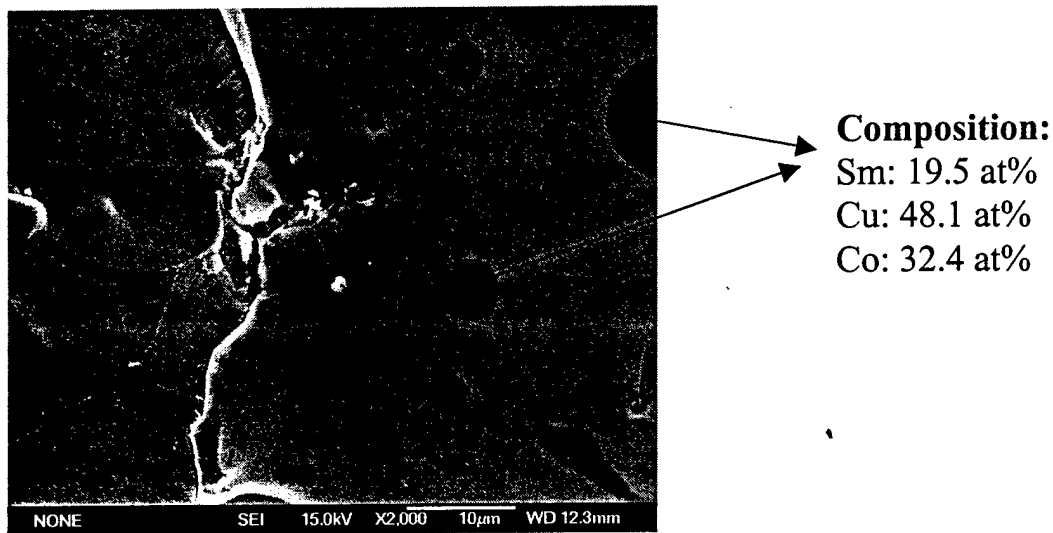


Fig. 1a. Typical SEM microstructure of cast SmCo_2Cu_3 magnet

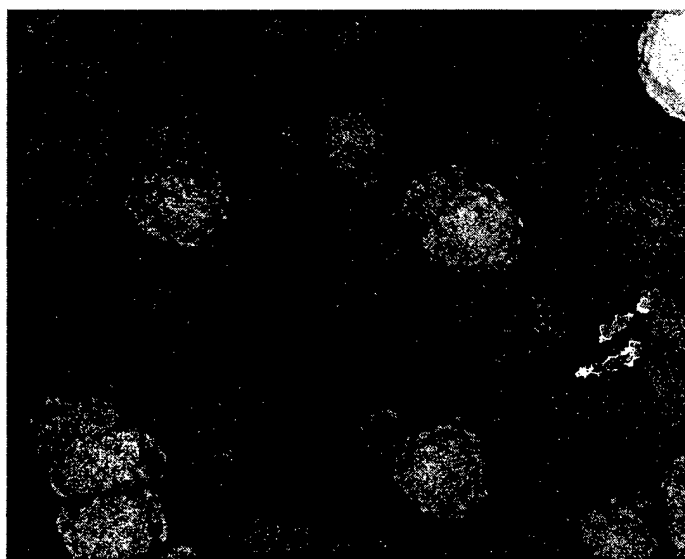


Fig. 1b. SEM image near grain boundary in cast SmCo_2Cu_3 magnet

Table 1 EDX chemical composition analysis of different particles in SmCo_2Cu_3 cast magnet

| Particle | Sm(at%) | Co (at%) | Cu (at%) | O (at%) |
|----------|---------|----------|----------|---------|
| 1 | 20.1 | 28.4 | 50.6 | 0.9 |
| 2 | 87 | 1.3 | 5.6 | 6.1 |
| 3 | 92 | 1.5 | 2.2 | 4.3 |
| 4 | 90.9 | 1.5 | 6.9 | 0.7 |
| 5 | 21.6 | 28.8 | 48 | 0.6 |
| 6 | 20.5 | 32.6 | 46.2 | 0.7 |

Since the installation of the high powered SEM, various types of methods in the production of nanoparticles can now be characterized. Conventional methods, such as deposition using standard sputtering techniques and chemical deposition are being investigated to produce nanoparticles. The SEM has been used to analyzed these various techniques, and improvements in the production of nanoparticles have been achieved. Below is a SEM image of TiO_2 film chemically deposited on a Si substrate.



TiO_2 clusters showing particles between 10 to 50 nm

Numerous programs in other departments including Material Science, Mechanical Engineering, Chemical Engineering, Geology, and Chemistry, have benefited from the presence of a modern SEM. In addition, the microscope has offered the type of educational training in state of the art microscopy and microanalysis skill needed for our student and postdoctoral fellows.

Upon completion of the installation, work proceeded in making the new Scanning Electron Microscope available for remote operation via the Internet. This has been accomplished and demonstrations were given to local High School Science teachers. Remote microscopy is an emerging technology that will make collaboration with distant colleges easier. The JEOL JSM 6335F is capable of remote operation anywhere an Ethernet connection is available.

FINAL REPORT
(DURIP 00) Scanning Electron Microscope (SEM)

AFOSR Award# F49620-00-1-0261
UD Account# 3-3-21-2570-03

April 1, 2000 – March 31, 2001


George C. Hadjipanayis
Principal Investigator

Submitted to:
Dr. Harold Weinstock
AFOSR/NE
801 N. Randolph Street
Room 732
Arlington, VA 22203-1977
(703) 696-8572


The purchase of a high power field emission Scanning Electron Microscope under DURIP FY2000 grant F49620-00-1-0261 has enhanced the Magnetism lab microanalysis capabilities. The Department of Physics at the University of Delaware purchased a state-of-the-art JEOL JSM 6335F Scanning Electron Microscope in the fall of 2000. The purchase price of \$235,000.00 was negotiated to include a high resolution Backscatter detector, motorized X and Y stage control, and infrared specimen chamber scope. The installation was completed in a timely manner and the instrument was turned over to the Magnetism Lab in November 2000.

The increased analytical capabilities of the Magnetism Lab, since the purchase of the field emission SEM, has made the observation of nanostructures in hard and soft magnetic material routine. Having a high power field emission SEM means the laborious sample preparation and extensive time required for TEM imaging is no longer required. Typically, the resolution of the SEM (1.5nm at 15KV) and its depth of field and microanalysis capabilities provide detailed information of samples that previously required sending them to various other facilities. One of the unique features of the field emission SEM is the ability to view elemental maps using EDS techniques. The relatively small beam size provides the resolution to map microstructures and to identify inhomogeneities in samples in the range of a few nanometers. Since the installation of the JSM 6335F, our research activities have increased. Quick screening of samples and their associated microstructure can now be performed in a timely manner.

Actual JSM6335F images

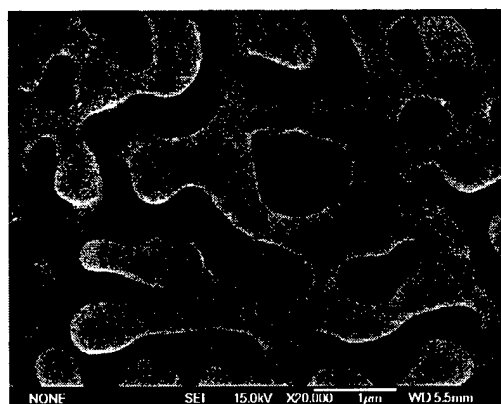

Porous Gold Microstructure

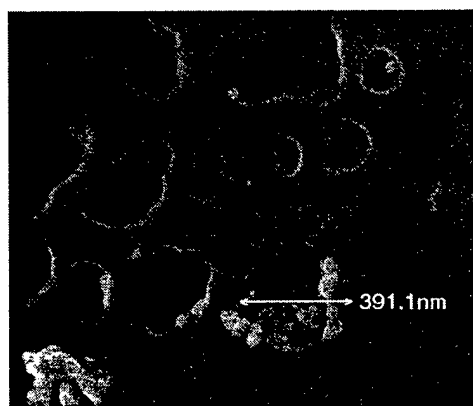

Silver Nanoparticles in
Polystyrene

One of the more recent undertakings of the UD magnetics lab has been the production of nanoparticles in a matrix using a novel sputtering based deposition technique. Shown below are SEM images of Fe nanoparticles in a BN matrix. Images at 30k magnification give particle sizes in the range of 20~100nm. At higher magnification (330k), it was possible to investigate clustering of nanoparticles. The image clearly shows the separation of the nanoparticles of a cluster enveloped by a possible passivation layer.

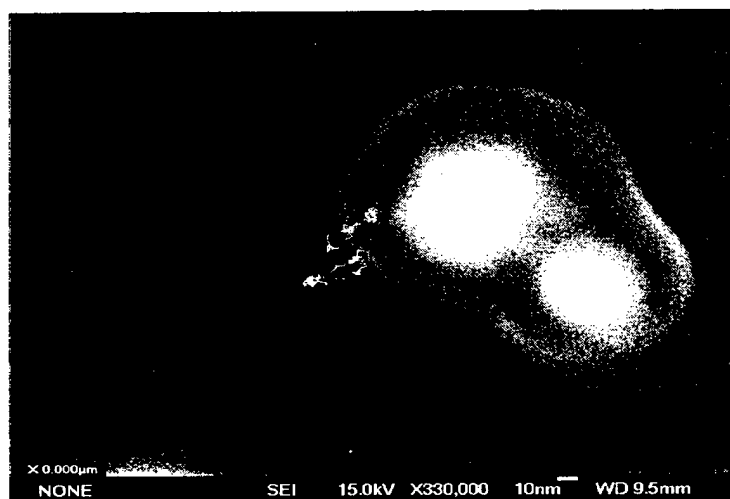
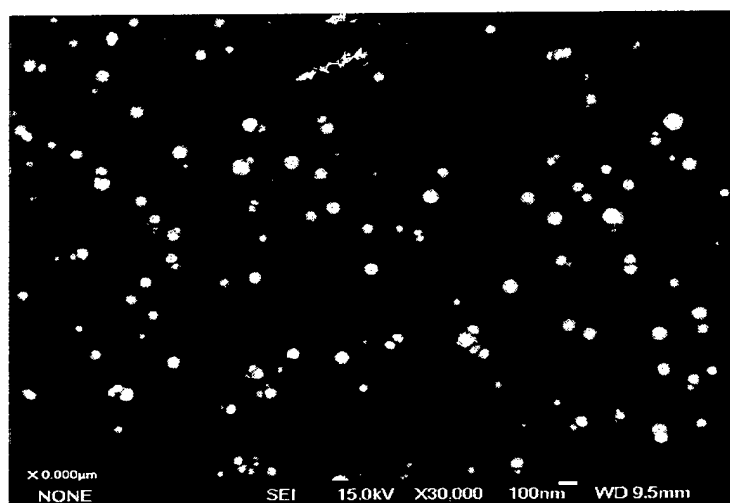


Figure 1 shows a typical SEM microstructure of cast SmCo_2Cu_3 magnets with EDS chemical composition analysis. There are many black dots in the right grain in figure 1a. EDS results found these dots had the same composition as matrix-1:5 phase, which means that these dots are originally from matrix phase. In higher magnification SEM images, some irregular particles were observed near the grain boundary (figure 1b). From EDS chemical composition analysis (table1), particle 1, particle 5, and particle 6 are the pieces of matrix phase, whereas particle 2, particle 3, and particle 4 show very high Sm content and considerable O content, which indicates that these particles are Sm oxides.

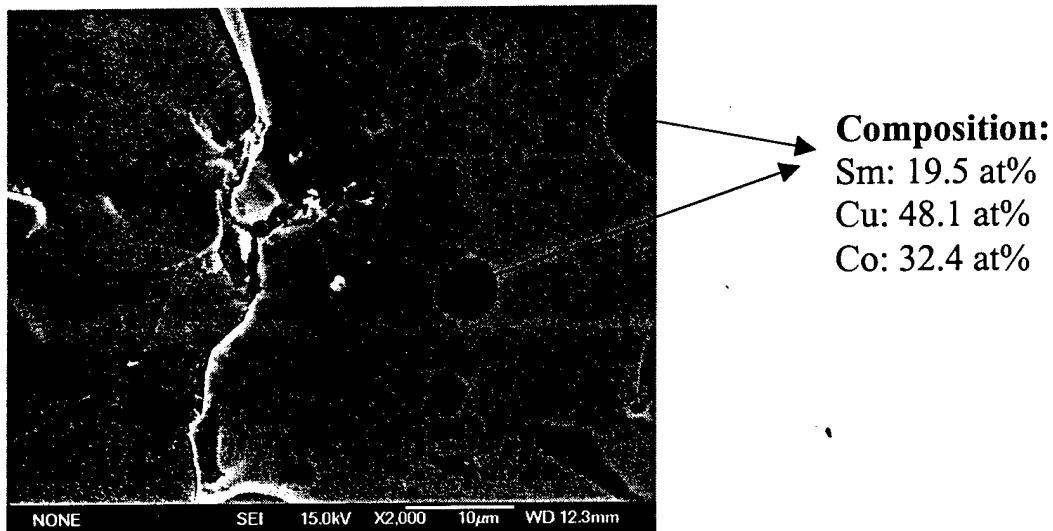


Fig. 1a. Typical SEM microstructure of cast SmCo_2Cu_3 magnet

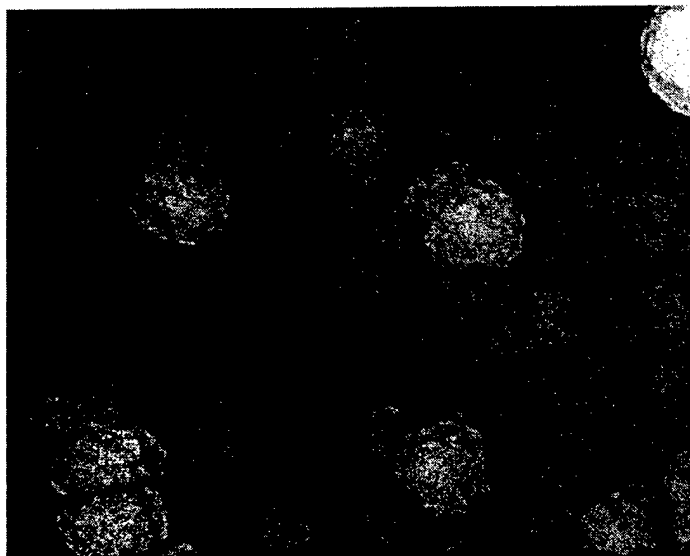


Fig. 1b. SEM image near grain boundary in cast SmCo_2Cu_3 magnet

Table 1 EDX chemical composition analysis of different particles in SmCo_2Cu_3 cast magnet

| Particle | Sm(at%) | Co (at%) | Cu (at%) | O (at%) |
|----------|---------|----------|----------|---------|
| 1 | 20.1 | 28.4 | 50.6 | 0.9 |
| 2 | 87 | 1.3 | 5.6 | 6.1 |
| 3 | 92 | 1.5 | 2.2 | 4.3 |
| 4 | 90.9 | 1.5 | 6.9 | 0.7 |
| 5 | 21.6 | 28.8 | 48 | 0.6 |
| 6 | 20.5 | 32.6 | 46.2 | 0.7 |

Since the installation of the high powered SEM, various types of methods in the production of nanoparticles can now be characterized. Conventional methods, such as deposition using standard sputtering techniques and chemical deposition are being investigated to produce nanoparticles. The SEM has been used to analyzed these various techniques, and improvements in the production of nanoparticles have been achieved. Below is a SEM image of TiO_2 film chemically deposited on a Si substrate.



TiO_2 clusters showing particles between 10 to 50 nm

Numerous programs in other departments including Material Science, Mechanical Engineering, Chemical Engineering, Geology, and Chemistry, have benefited from the presence of a modern SEM. In addition, the microscope has offered the type of educational training in state of the art microscopy and microanalysis skill needed for our student and postdoctoral fellows.

Upon completion of the installation, work proceeded in making the new Scanning Electron Microscope available for remote operation via the Internet. This has been accomplished and demonstrations were given to local High School Science teachers. Remote microscopy is an emerging technology that will make collaboration with distant colleges easier. The JEOL JSM 6335F is capable of remote operation anywhere an Ethernet connection is available.

REPORT DOCUMENTATION PAGE

AFRL-SR-BL-TR-01-

Public reporting burden for this collection of information is estimated to average 1 hour per response, including gathering and maintaining the data needed, and completing and reviewing the collection of information. Send collection of information, including suggestions for reducing this burden, to Washington Headquarters Service, Paperwork Project, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Project, Suite 1204, Arlington, VA 22202-4302.

if this
erson

| | | | |
|--|---|--|---|
| 1. AGENCY USE ONLY (Leave blank) | | 2. REPORT DATE | 3. REPORT TYPE AND DATES COVERED 01 Apr 00 to 31 Mar 01 Final |
| 4. TITLE AND SUBTITLE (DURIP 00) Scanning Electron Microscope (SEM) | | | 5. FUNDING NUMBERS 61103D 3484/US |
| 6. AUTHOR(S) Professor Hadjipansyis | | | |
| 7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) Univ of Delaware 210 Hullihen Hall Newark, DE 19716-1551 | | | 8. PERFORMING ORGANIZATION REPORT NUMBER |
| 9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES) AFOSR/NE 801 North Randolph Street Rm 732 Arlington, VA 22203-1977 | | | 10. SPONSORING/MONITORING AGENCY REPORT NUMBER F49620-00-1-0261 |
| 11. SUPPLEMENTARY NOTES | | | |
| 12a. DISTRIBUTION AVAILABILITY STATEMENT APPROVAL FOR PUBLIC RELEASE; DISTRIBUTION UNLIMITED | | | 12b. DISTRIBUTION CODE 130-12. DISTRIBUTION IS UNLIMITED. |
| 13. ABSTRACT (Maximum 200 words) Upon completion of the installation, work proceeded in making the new Scanning Electron Microscope available for remote operation via the Internet. This has been accomplished and demonstrations were given to local High School Science teachers. Remote microscopy is an emerging technology that will make collaboration with distant colleges easier. The JEOL JSM 6335F is capable of remote operation anywhere an Ethernet connection is available. | | | |
| 14. SUBJECT TERMS | | | 15. NUMBER OF PAGES |
| | | | 16. PRICE CODE |
| 17. SECURITY CLASSIFICATION OF REPORT UNCLASSIFIED | 18. SECURITY CLASSIFICATION OF THIS PAGE UNCLASSIFIED | 19. SECURITY CLASSIFICATION OF ABSTRACT UNCLASSIFIED | 20. LIMITATION OF ABSTRACT UL |